
# seca 384/385

## Instructions for Use

17-10-07-688-100a\_2026-09S
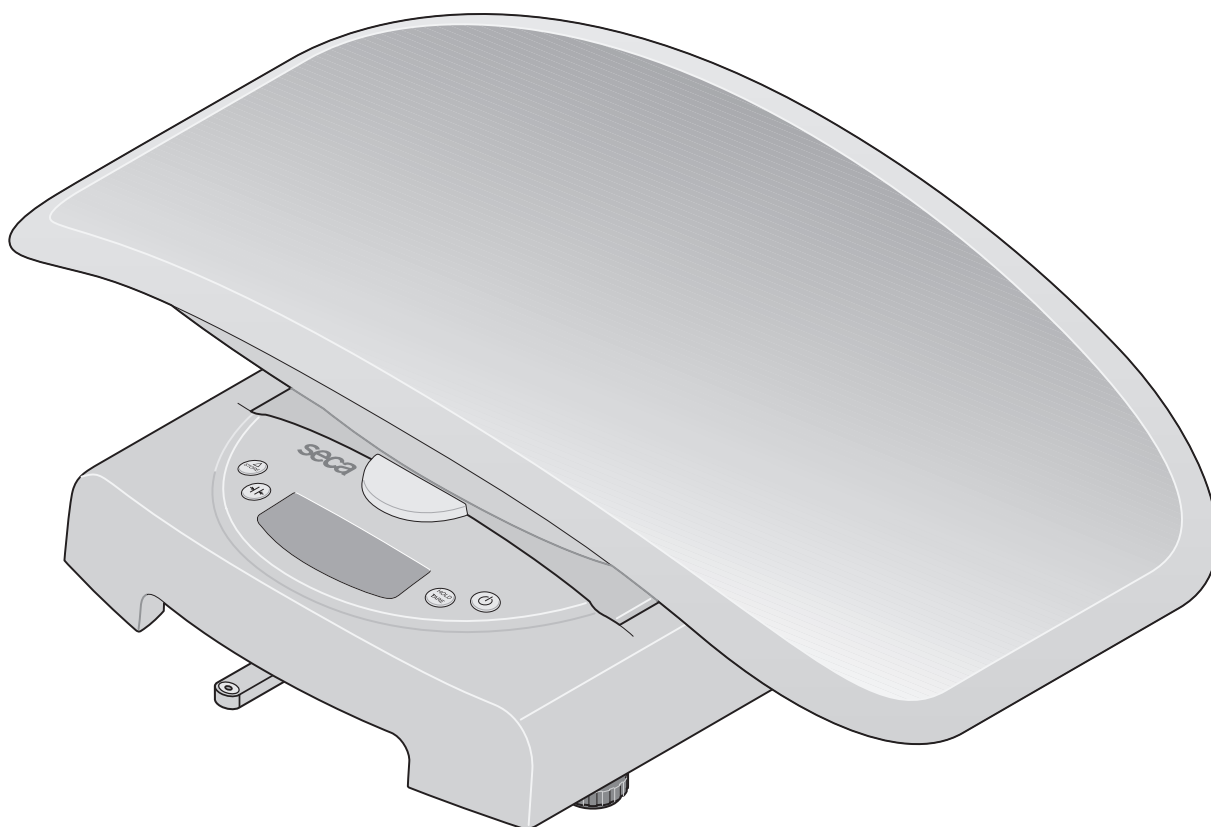
|                                                             |            |
|-------------------------------------------------------------|------------|
| <b>Deutsch – Gebrauchsanweisung .....</b>                   | <b>3</b>   |
| <b>English – Instructions for Use.....</b>                  | <b>24</b>  |
| <b>Français – Mode d’emploi.....</b>                        | <b>45</b>  |
| <b>Italiano – Istruzioni per l’uso .....</b>                | <b>66</b>  |
| <b>Español – Instrucciones de uso .....</b>                 | <b>87</b>  |
| <b>Português – Manual de instruções de utilização .....</b> | <b>108</b> |
| <b>Polski – Instrukcja obsługi.....</b>                     | <b>129</b> |
| <b>Dansk – Brugsanvisning.....</b>                          | <b>150</b> |
| <b>Svenska – Bruksanvisning .....</b>                       | <b>170</b> |
| <b>Norsk – Bruksanvisning.....</b>                          | <b>191</b> |
| <b>Suomi – Käyttöohje .....</b>                             | <b>212</b> |
| <b>Nederlands – Gebruiksaanwijzing .....</b>                | <b>232</b> |
| <b>Ελληνικά – Οδηγίες χρήσης.....</b>                       | <b>253</b> |
| <b>Česky – Návod k použití .....</b>                        | <b>275</b> |

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# CONTENTS

|                                                               |           |                                                                 |           |
|---------------------------------------------------------------|-----------|-----------------------------------------------------------------|-----------|
| <b>1 About this document .....</b>                            | <b>25</b> | <b>11 Technical data.....</b>                                   | <b>41</b> |
| 1.1 Representation in text.....                               | 25        | 11.1 General technical data.....                                | 41        |
| 1.2 Representation in diagrams .....                          | 25        | 11.2 Weight measurement .....                                   | 42        |
| <b>2 Description of device .....</b>                          | <b>25</b> | <b>12 Optional accessories and spare parts.....</b>             | <b>43</b> |
| 2.1 Intended use.....                                         | 25        | <b>13 Disposal .....</b>                                        | <b>43</b> |
| 2.2 Patient target group.....                                 | 25        | 13.1 Disposing of the device .....                              | 43        |
| 2.3 User qualification .....                                  | 25        | 13.2 Disposing of batteries and rechargeable<br>batteries ..... | 43        |
| 2.4 Contraindications.....                                    | 26        | <b>14 Warranty .....</b>                                        | <b>44</b> |
| 2.5 Clinical benefit .....                                    | 26        | <b>15 Declaration of conformity .....</b>                       | <b>44</b> |
| 2.6 Functional description.....                               | 26        |                                                                 |           |
| Measuring weight .....                                        | 26        |                                                                 |           |
| Functional description.....                                   | 26        |                                                                 |           |
| <b>3 Safety precautions.....</b>                              | <b>26</b> |                                                                 |           |
| 3.1 Safety information in these instructions for<br>use ..... | 26        |                                                                 |           |
| 3.2 Basic safety information.....                             | 26        |                                                                 |           |
| Handling device .....                                         | 26        |                                                                 |           |
| Preventing injuries and infections .....                      | 27        |                                                                 |           |
| Preventing damage to device.....                              | 28        |                                                                 |           |
| Handling measuring results.....                               | 28        |                                                                 |           |
| Handling packaging material.....                              | 29        |                                                                 |           |
| Handling batteries and rechargeable<br>batteries .....        | 29        |                                                                 |           |
| <b>4 Overview .....</b>                                       | <b>30</b> |                                                                 |           |
| 4.1 Scale .....                                               | 30        |                                                                 |           |
| 4.2 Display.....                                              | 31        |                                                                 |           |
| 4.3 Labels.....                                               | 31        |                                                                 |           |
| <b>5 Starting up device.....</b>                              | <b>33</b> |                                                                 |           |
| 5.1 Scope of delivery .....                                   | 33        |                                                                 |           |
| 5.2 Transporting the device .....                             | 34        |                                                                 |           |
| 5.3 Setting up and aligning the device .....                  | 34        |                                                                 |           |
| 5.4 Assembling/disassembling the tray .....                   | 34        |                                                                 |           |
| 5.5 Inserting batteries .....                                 | 35        |                                                                 |           |
| <b>6 Operation.....</b>                                       | <b>35</b> |                                                                 |           |
| 6.1 Switching the device on and off .....                     | 35        |                                                                 |           |
| Switching on the device.....                                  | 35        |                                                                 |           |
| Switching off the device .....                                | 35        |                                                                 |           |
| 6.2 Measuring weight .....                                    | 35        |                                                                 |           |
| Weighing .....                                                | 36        |                                                                 |           |
| Taring additional weight (Tare) .....                         | 36        |                                                                 |           |
| Displaying measuring result permanently<br>(Hold).....        | 37        |                                                                 |           |
| Using the differential value function .....                   | 37        |                                                                 |           |
| Switching the weighing range .....                            | 38        |                                                                 |           |
| <b>7 Hygiene treatment .....</b>                              | <b>39</b> |                                                                 |           |
| 7.1 Cleaning .....                                            | 39        |                                                                 |           |
| 7.2 Disinfecting.....                                         | 39        |                                                                 |           |
| 7.3 Sterilizing .....                                         | 39        |                                                                 |           |
| <b>8 Function check .....</b>                                 | <b>40</b> |                                                                 |           |
| <b>9 Troubleshooting.....</b>                                 | <b>40</b> |                                                                 |           |
| 9.1 General faults .....                                      | 40        |                                                                 |           |
| <b>10 Servicing.....</b>                                      | <b>41</b> |                                                                 |           |

# 1 ABOUT THIS DOCUMENT

## 1.1 Representation in text

| Symbol   | Description                                |
|----------|--------------------------------------------|
| ✓        | Requirement for actions                    |
| ►        | Action                                     |
| 1.<br>2. | Actions with specified sequence            |
| a)<br>b) | Steps of an action with specified sequence |
| ⇒        | Result of an action                        |
| •<br>•   | First level of a list                      |
| –<br>–   | Second level of a list                     |

## 1.2 Representation in diagrams

| Symbol | Description                                                                                                |
|--------|------------------------------------------------------------------------------------------------------------|
| ➡      | Indicates directions of movement of the device or of device components                                     |
| ✓      | <ul style="list-style-type: none"> <li>• Correct action</li> <li>• Correct result of action</li> </ul>     |
| ✗      | <ul style="list-style-type: none"> <li>• Incorrect action</li> <li>• Incorrect result of action</li> </ul> |

# 2 DESCRIPTION OF DEVICE

## 2.1 Intended use

The electronic baby and child scale supports physicians in decision-making regarding diagnosis or therapy, based on the patient's weight or weight changes.

## 2.2 Patient target group

The device is intended for babies, toddlers and children who do not exceed the maximum capacity of the scale.

## 2.3 User qualification

The device may only be operated by persons with formal training in the healthcare sector or in medicine.

## 2.4 Contraindications

No contraindications are known.

## 2.5 Clinical benefit

The electronic baby and child scale supports physicians in decision-making regarding diagnosis or therapy, based on the measured weight. With that, the electronic baby and child scale only provides an indirect clinical benefit. The physician must perform additional examinations for an accurate diagnosis.

## 2.6 Functional description

Weight calculation is performed by four load cells.

The tray of the scale can be taken off to weigh children who can stand independently.

# 3 SAFETY PRECAUTIONS

## 3.1 Safety information in these instructions for use

### **DANGER!**

Used to identify an extremely hazardous situation. If you fail to take note of this information, serious irreversible or fatal injuries will occur.

### **WARNING!**

Used to identify an extremely hazardous situation. If you fail to take note of this information, serious irreversible or fatal injuries may result.

### **CAUTION!**

Used to identify a hazardous situation. If you fail to take note of this information, minor to moderate injuries may result.

### **NOTICE!**

Used to identify possible incorrect usage of the device. If you fail to take note of this information, you may damage the device, or the measuring results may be incorrect.

### **NOTE**

Includes additional information about use of the device.

## 3.2 Basic safety information

### **Handling device**

- ▶ Please take note of the information in these instructions for use.
- ▶ Keep the instructions for use in a safe place. The instructions for use are a component of the device and must be available at all times.
- ▶ In the interest of patient safety, you and your patients are obliged to report serious events that occur in connection with this product to the manufacturer and to the authority responsible in your country.

**DANGER!****Risk of explosion**

- ▶ Do not use the device in an environment in which one of the following gases has accumulated:
  - oxygen
  - flammable anesthetics
  - other flammable substances/air mixtures

**CAUTION!****Patient hazard, damage to device**

- ▶ Additional devices connected to medical electrical devices must provide evidence of compliance with the relevant IEC or ISO standards (e.g. IEC 60950 for data-processing devices). Furthermore, all configurations must comply with the requirements of standards for medical systems (see IEC 60601-1-1 or Section 16 of edition 3.1 of IEC 60601-1 respectively). Anyone connecting additional devices to medical electrical devices is considered a system configurator and therefore responsible for ensuring that the system complies with the requirements of standards for systems. This also applies to additional devices recommended by seca. Your attention is drawn to the fact that local laws take precedence over the above-mentioned requirements of standards. In the event of any queries, please contact your local specialist dealer or Technical Service.

**CAUTION!****Patient hazard, damage to device**

- ▶ Have servicing carried out regularly as described in the relevant section of this document.
- ▶ Technical modifications may not be made to the device. The device does not contain any user-serviceable parts. Only have servicing and repairs performed by an authorized seca service partner. You can find your local service partner at [www.seca.com](http://www.seca.com).
- ▶ Use only seca original accessories and spare parts, otherwise seca will not grant any warranty.

**CAUTION!****Patient hazard, malfunction**

- ▶ Keep medical electrical devices such as high-frequency surgical equipment at a minimum distance of approx. 1 meter to prevent faulty measurements.
- ▶ Keep high-frequency equipment such as cellphones at a minimum distance of approx. 1 meter to prevent faulty measurements.
- ▶ The actual transmission output of HF devices may require minimum distances of more than 1 meter. For details, go to [www.seca.com](http://www.seca.com).

**Preventing injuries and infections****WARNING!****Injury from falling**

Baby scales are usually placed on elevated work surfaces. If the baby were to fall from the work surface, it could suffer serious irreversible or fatal injuries.

- ▶ Route the connector cables (if present) to prevent users from tripping over them.
- ▶ Ensure that the device is steady and level.
- ▶ Never leave a baby unsupervised.

## Preventing damage to device

### **WARNING!**

#### **Risk of infection**

- ▶ Before and after every measurement, wash your hands to reduce the risk of cross-contamination and nosocomial infections.
- ▶ Should the patient have contagious diseases, subject the device to a hygiene treatment as described in the relevant section of these instructions for use immediately afterwards.
- ▶ Ensure that the patient does not have any open wounds or infectious skin alterations which may come into contact with the device.
- ▶ Subject the device to a hygiene treatment at regular intervals.

### **NOTICE!**

#### **Damage to device**

- ▶ Ensure that fluids never get inside the device. These can destroy the electronics.
- ▶ For devices with power supply operation: Switch off the device before you disconnect the power supply unit from the power supply socket.
- ▶ For devices with power supply operation: If the device is not to be used for an extended period, disconnect the power supply unit from the power supply socket. Only then is the device de-energized.
- ▶ For devices with battery or rechargeable battery operation: If you are not using the device for an extended period of time, remove batteries or battery packs. Only then is the device de-energized.
- ▶ Do not drop the device.
- ▶ Do not subject the device to shocks or vibrations.
- ▶ Perform a function check before each use as described in the corresponding section in this document. Do not operate the device if it is not working properly or is damaged.
- ▶ Do not place the device in direct sunlight and ensure that it is not placed in the direct proximity of a heat source. The excessive temperatures could damage the electronics.
- ▶ Avoid rapid temperature fluctuations. If the device is transported so that a temperature difference of over 20 °C occurs, the device must be left to stand for at least 2 hours before it is switched on, otherwise condensation may form; this may damage the electronics.
- ▶ Use the device only in the intended ambient conditions outlined in the section entitled "Technical data".
- ▶ Store the device only in the intended storage conditions outlined in the section entitled "Technical data".
- ▶ Use only cleaning agents and disinfectants which match the details in the section entitled "Hygiene treatment".

### **NOTICE!**

#### **Damage to device**

- ▶ Ensure that maximum capacity is not exceeded.

## Handling measuring results

### **NOTICE!**

#### **Inconsistent measuring results**

- ▶ Before you save and continue using measured values determined using this device (e.g. in seca software or in an EMR system), ensure that the measured values are plausible.
- ▶ If measured values have been sent to seca software or to an EMR system, ensure before continuing to use them that the measured values are plausible and assigned to the correct patient.

### Handling packaging material

#### **WARNING!**

##### **Risk of suffocation**

Packaging materials made of plastic film (e.g. bags) present a risk of suffocation.

- ▶ Store packaging material out of the reach of children.
- ▶ If the original packaging material is no longer available, only use plastic bags with safety holes to reduce the risk of suffocation. Use recyclable materials if possible.

#### **NOTE**

Store the original packaging material for future use (e.g. returning for servicing).

### Handling batteries and rechargeable batteries

#### **WARNING!**

##### **Injury**

Batteries and rechargeable batteries contain harmful substances which may explode and so escape if not handled properly.

- ▶ Do not try to recharge batteries.
- ▶ Do not expose batteries/rechargeable batteries to heat.
- ▶ Do not burn batteries/rechargeable batteries.
- ▶ If harmful substances are leaking out, avoid contact with the skin, eyes, and mucous membranes. Rinse affected areas of the body with plenty of clean water and seek medical help at once.


#### **NOTICE!**

##### **Damage to device and malfunction due to incorrect handling**

- ▶ Use only the type of battery/rechargeable battery quoted in this document.
- ▶ Always replace all batteries/rechargeable batteries at the same time.
- ▶ Do not short-circuit batteries/rechargeable batteries.
- ▶ If you do not use the device for an extended period of time, remove the batteries/rechargeable batteries. This prevents acid from leaking into the device.
- ▶ If harmful substances have penetrated the device, do not continue using it. Have the device checked by an authorized service partner and repaired if necessary.


## 4 OVERVIEW

### 4.1 Scale
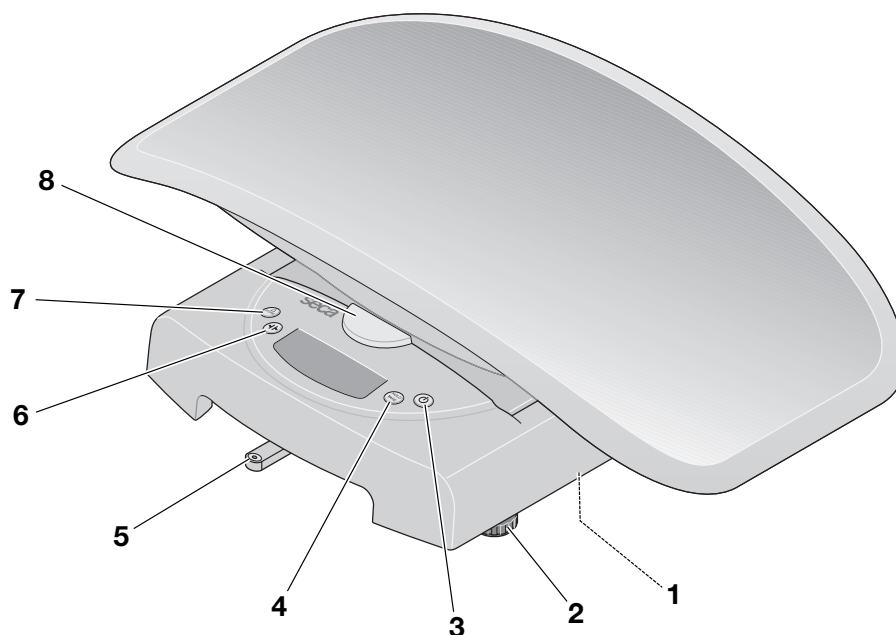
| Item | Element                                                                             | Function                                                                                                                                              |
|------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------|
| 1    | Battery compartment                                                                 | Takes batteries (6 x type AA/1.5 V)                                                                                                                   |
| 2    | Foot screw                                                                          | Used for precise alignment of the scale (4 pcs.)                                                                                                      |
| 3    |  | Used for switching the device on and off                                                                                                              |
| 4    |  | <ul style="list-style-type: none"> <li>Press briefly: Activate <b>Hold</b> function</li> <li>Press and hold: Activate <b>Tare</b> function</li> </ul> |
| 5    | Spirit level                                                                        | Indicates whether the device is horizontal                                                                                                            |
| 6    |  | Used for switching weighing range                                                                                                                     |
| 7    |  | Saves the weight value (prior to food intake → <a href="#">Using the differential value function, page 37</a> )                                       |
| 8    | Release button                                                                      | Used to release the tray                                                                                                                              |

## 4.2 Display
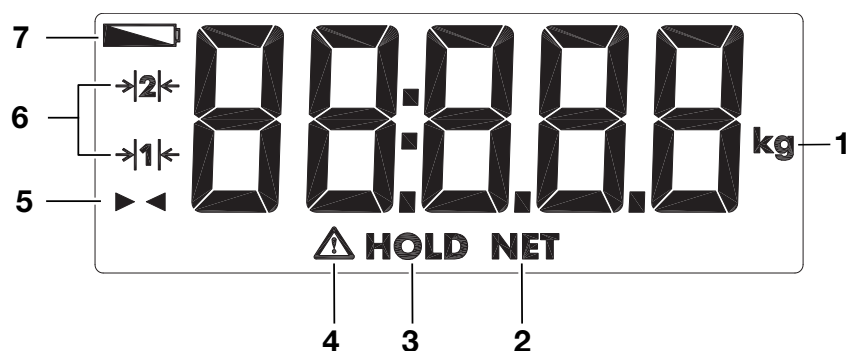
| Item | Symbol       | Meaning                            |
|------|--------------|------------------------------------|
| 1    | kg           | Unit in which weight is displayed  |
| 2    | NET          | <b>Tare</b> function active        |
| 3    | HOLD         | <b>Hold</b> function active        |
| 4    |              | Non-verifiable function active     |
| 5    |              | Differential value function active |
| 6    | →1 ←<br>→2 ← | Weighing range currently in use    |
| 7    |              | Batteries are low                  |

## 4.3 Labels

| Markings on the device and on the type plate |                                                                                                                                                                             |
|----------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Text/symbol                                  | Meaning                                                                                                                                                                     |
|                                              | Name and address of manufacturer, date of manufacture                                                                                                                       |
| UDI                                          | Unique Device Identification<br>(product identification number as per Directive (EU) 2017/745)                                                                              |
| REF                                          | Article number                                                                                                                                                              |
| SN                                           | Serial number                                                                                                                                                               |
| GAL                                          | Value in m/s <sup>2</sup> (model-dependent) <ul style="list-style-type: none"> <li>Gravitational acceleration on earth</li> <li>Depends on the intended location</li> </ul> |
| ProdID                                       | seca product identification number                                                                                                                                          |
| Approval Type                                | Type designation of design approval                                                                                                                                         |
| Mat.No.                                      | Variant number                                                                                                                                                              |
|                                              | Follow instructions for use                                                                                                                                                 |


| Markings on the device and on the type plate                                                                                                                               |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Text/symbol                                                                                                                                                                | Meaning                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |
|                                                                                           | Medical electrical device, Type B                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                            |
| IP20                                                                                                                                                                       | Type of protection in accordance with IEC 60529: <ul style="list-style-type: none"> <li>• Protection against ingress of solid foreign objects with a diameter of over 12.5 mm</li> <li>• Protection against access with fingers</li> <li>• No protection from ingress of water</li> </ul>                                                                                                                                                                                                                                                                    |
| e                                                                                                                                                                          | Value in units of mass<br>Used to classify and verify a scale                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                |
|                                                                                           | Class III scale as per directive 2014/31/EU                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |
|  M24 <sup>0102</sup> <sub>0123</sub>                                                      | Device complies with EU directives <ul style="list-style-type: none"> <li>• <b>M</b>: Conformity label in compliance with directive 2014/31/EU for non-automatic scales (verified models)</li> <li>• <b>26</b>: (Example: 2026) Year in which declaration of conformity was implemented and the CE-label was applied (verified models)</li> <li>• <b>0102</b>: Notified Body for Metrology (verified models)</li> <li>• <b>0123</b>: Notified Body for Medical Devices</li> </ul>                                                                            |
|                                                                                           | Medical device in accordance with Regulation (EU) 2017/745                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                   |
|  M24 <sup>xxxx</sup> <sub>yyyy</sub>                                                    | Device complies with United Kingdom directives <ul style="list-style-type: none"> <li>• <b>M</b>: Conformity label in compliance with UK directive SI 2016 no. 1152 for non-automatic scales (NAWIR) (verified models)</li> <li>• <b>24</b>: (Example: 2024) Year in which declaration of conformity was implemented and the UKCA label was applied (verified models)</li> <li>• <b>xxxx</b>: Notified Body for Medical Devices of the United Kingdom</li> <li>• <b>yyyy</b>: Notified Body for Metrology of the United Kingdom (verified models)</li> </ul> |
| <br> | Importer/representative in the United Kingdom:<br>seca Ltd<br>40 Barn Street<br>B5 5QB Birmingham<br>United Kingdom                                                                                                                                                                                                                                                                                                                                                                                                                                          |
| <br> | Importer/representative in Switzerland:<br>seca ag (schweiz)<br>Medizinische Waagen und Messsysteme<br>Schön matt Str. 2<br>CH-4153 REINACH                                                                                                                                                                                                                                                                                                                                                                                                                  |
|                                                                                         | Do not dispose of device in household waste                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |

| Markings on the packaging                                                           |                                                                              |
|-------------------------------------------------------------------------------------|------------------------------------------------------------------------------|
| Symbol                                                                              | Meaning                                                                      |
|  | Protect from moisture                                                        |
|  | Arrows indicate top of product<br>Transport and store in an upright position |
|  | Fragile, do not throw or drop                                                |

| Markings on the packaging                                                         |                                                                  |
|-----------------------------------------------------------------------------------|------------------------------------------------------------------|
| Symbol                                                                            | Meaning                                                          |
|  | Permitted min. and max. temperature for transport and storage    |
|  | Permitted min. and max. humidity for transport and storage       |
|  | Permitted min. and max. air pressure for transport and storage   |
|  | Open packaging here                                              |
|  | Packaging material can be disposed of through recycling programs |

5     **STARTING UP DEVICE**


5.1    **Scope of delivery**
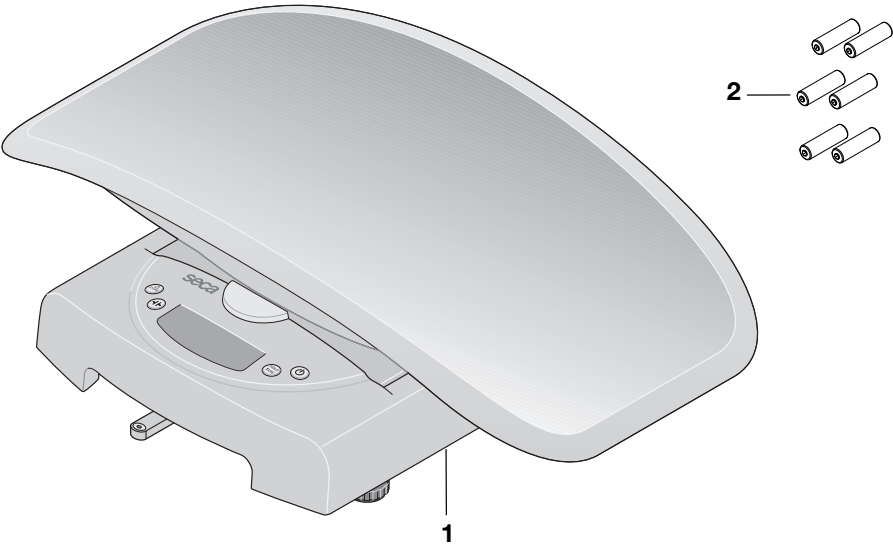
| Item      | Component                 | Pcs. |
|-----------|---------------------------|------|
| 1         | Baby scale                | 1    |
| 2         | Batteries (type AA/1.5 V) | 6    |
| Not shown | Instructions for use      | 1    |

17-10-07-688-100a\_2026-09S

## 5.2 Transporting the device

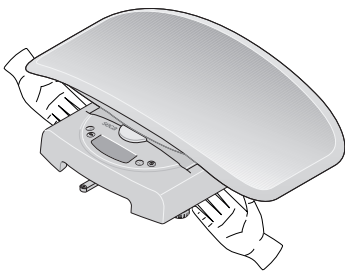
### NOTICE!

#### Damage to device from incorrect handling

The tray is not intended for transport of the scale. The tray can become deformed or come off the scale as a result of high strain.

- When transporting the scale, hold it only by the housing.
- Lift the scale by its housing as shown in the illustration.

## 5.3 Setting up and aligning the device
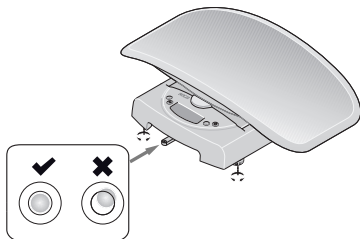
### CAUTION!

#### Faulty measurement as a result of devices being set up incorrectly

- Place the device on a stable, even work surface or on a seca cart (Baby Scale Cart).
- Align the device before using it and every time it changes location.
- Align the device only without a load on it.


### WARNING!

#### Faulty measurement as a result of force shunt


- Make sure that the scale does not touch the wall or other objects during measurement.
1. Place the scale on a firm, level surface.
  2. Align the scale by turning the foot screws. The air bubble in the spirit level must be located in the center of the circle.
  3. Ensure that only the feet of the scale are in contact with the ground.

## 5.4 Assembling/disassembling the tray
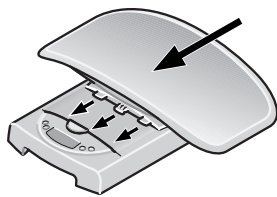
### Assembling the tray


1. Place the tray as shown on the weighing platform.
  2. Keep pulling the tray toward the front until the tray audibly engages and the release button is in the upper position.
- ⇒ The tray is firmly locked to the housing.
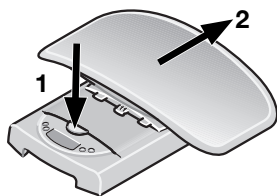
### Disassembling the tray

1. Firmly press the release button and push the tray to the back.

## 5.5 Inserting batteries

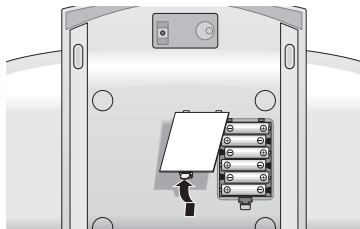✓ Six batteries (type AA/1.5 V)

1. Take off the tray (→ [Assembling/disassembling the tray, page 34](#)).
2. Turn the scale carefully to ensure that the underside is accessible.
3. Open the battery compartment.

### NOTICE!

#### Damage to device

The device may overheat and become damaged.

- Ensure the correct polarity of the batteries (markings in the battery compartment).

4. Place the batteries in the battery compartment.
5. Close the battery compartment.
6. Carefully turn the scale over.
7. Assemble the tray (→ [Assembling/disassembling the tray, page 34](#)).

## 6 OPERATION

### WARNING!

#### Personal injury, faulty measurement

- Before using the device each time, perform a function check as described in the corresponding section of these instructions for use.

### 6.1 Switching the device on and off

#### Switching on the device

✓ Ensure that the device has no load.

- Press the  key.
  - ⇒ **SECA** is displayed.
  - ⇒ All the elements of the display are shown briefly.
  - ⇒ The device is ready for operation when **0.000** is displayed.

#### Switching off the device

- Press the ON/OFF button .

### NOTE

After it is switched on, the scale remains ready for weighing for approx. 20 seconds. If a weighing operation was performed, the display value remains for approx. 2 minutes until the scale switches off automatically.

### 6.2 Measuring weight

### WARNING!

#### Injury from falling

Baby scales are usually placed on elevated work surfaces. If the baby were to fall from the work surface, it could suffer serious irreversible or fatal injuries.

- Ensure that the device is steady and level.
- Never leave a baby unsupervised.

## Weighing

1. Select the weighing range (→ [Switching the weighing range, page 38](#)).

### Weighing babies


- ✓ The device has no load.

1. Press the  key.
2. Place the baby on the scale.
3. Make sure that the child is lying as still as possible.  
⇒ The display flashes until a stable weight is measured.
4. Read off the measuring result.

### Weighing toddlers and children

- ✓ The infant, toddler or child can stay still without assistance.

1. Take off the tray (→ [Assembling/disassembling the tray, page 34](#)).
2. Make sure that the weighing platform is placed on a stable surface.
- ✓ The device has no load.
3. Press the  key.
4. Ask the child to step onto the scale.
5. Make sure that the child is standing still.  
⇒ The display flashes until a stable weight is measured.
6. Read off the measuring result.
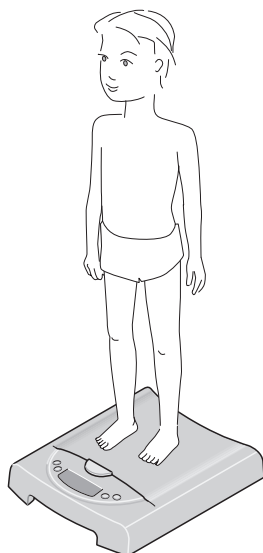
## Taring additional weight (Tare)

Using the **Tare** function, you can prevent any additional weight (e.g. a towel or a cushioned mat on the weighing surface) from affecting the measuring result.

### WARNING!

#### Faulty measurement with the Tare function

If the **Tare** function is used incorrectly, the weight of the patient may not be determined correctly.

- Perform the steps as described to prevent faulty measurements.

### WARNING!

#### Faulty measurement

If an additional weight is touching the surface on which the device is standing, the weight will not be measured correctly.

- Make sure that additional weights are positioned on the tray only or (if the child is standing) on the weighing platform without any overhang.

### NOTE

The maximum weight which can be displayed is reduced by the weight of the objects already placed on the scale.

- ✓ The device has no load.

1. Press the **START** key.
2. Place the additional weight on the scale.
3. Press and hold the **TARE** key.
4. Wait until the display stops flashing and **0.000** is displayed instead.
5. Place the baby on the scale.

6. Make sure that the child is lying as still as possible.  
⇒ The display flashes until a stable weight is measured.
7. Read off the measuring result.  
⇒ The additional weight is deducted automatically.

**NOTE**

Now, you can perform as many measurements as you like with this additional weight.

8. To deactivate the **Tare** function, press and hold the **TARE** key.

**NOTE**

When you switch the device off, the **Tare** function is deactivated automatically.

### Displaying measuring result permanently (Hold)

When you activate the **Hold** function, the weight value continues to be displayed after the load is removed from the scale. This allows you to take care of the baby before recording the weight value.

- ✓ The device has no load.

1. Press the key.
2. Place the baby on the scale.
3. Make sure that the child is lying as still as possible.  
⇒ The display flashes until a stable weight is measured.
4. Briefly press the **HOLD/TARE** key.  
⇒ The display flashes until a stable weight is measured.  
⇒ The weight is displayed permanently.  
⇒ The message **HOLD** is displayed.
5. Read off the measuring result.
6. To deactivate the **Hold** function, briefly press the **HOLD/TARE** key.

**NOTE**

The **Hold** function can also be used if the **Tare** function had been performed previously

**NOTE**

The value is displayed until the scale switches off automatically.

### Using the differential value function

You can use the differential value function to save a single weight value and perform differential calculations. For example, you can record how much milk the baby ingests during nursing or display the baby's weight gain/increase since its last weighing.

**Prior to food intake**

- ✓ The device has no load.

1. Press the key.
2. Place the baby on the scale.
3. Make sure that the child is lying as still as possible.  
⇒ The display flashes until a stable weight is measured.
4. Keep the **STORE** key depressed for approx. 2 seconds.  
⇒ The message **Sto** is displayed.  
⇒ The weight value is saved and retained after the scale is switched off.

**After food intake**

- ✓ The device has no load.

1. Press the key.
2. Place the baby on the scale.
3. Make sure that the child is lying as still as possible.  
⇒ The display flashes until a stable weight is measured.

4. Briefly press the **STORE** key.  
 ⇒ The ► ◀ symbol is displayed together with the differential value (current weight minus stored weight).
5. To switch off the differential value function, press the **STORE** key.

**NOTE**  
 When you save a new weight value, the previously saved value will be overwritten.


**NOTE**  
 The value is displayed until the scale switches off automatically.

**Switching the weighing range**

The scale has two weighing ranges. In weighing range 1 (→|1|←), graduations in the weight display are finer. In weighing range 2 (→|2|←), you can use the maximum capacity of the scale.

| Weighing range | Maximum capacity                     | Graduation |
|----------------|--------------------------------------|------------|
| 1              | 10 kg                                | 5 g        |
| 2              | 20 kg (seca 384)<br>50 kg (seca 385) | 10 g       |

- ✓ The device has no load.
1. Press the  key.  
 ⇒ The scale is working within the last weighing range selected.
  2. To switch the weighing range, press the key .  
 ⇒ The →|1|← symbol or →|2|← in the display demarcates the active weighing range.  
 ⇒ The active weighing range remains active after switch-off.

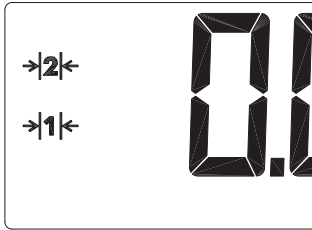**NOTE**  
 You can switch from range 1 to range 2 even if there is a load on the scale. Exception: A special function is activated (e.g. tare). For technical reasons, it is only possible to switch from range 2 to range 1 if there is a load on the scale. No special function may be activated here either.

## 7 HYGIENE TREATMENT

### **WARNING!** **Electric shock**

The device is not de-energized when the on/off button is pressed and the display goes out. Use of fluids on the device may cause an electric shock.

- ▶ Before each hygiene treatment, ensure that the device is switched off.
- ▶ Take the batteries/battery packs out of the device before every hygiene treatment.
- ▶ Ensure that no fluids penetrate the device.

### **WARNING!** **Risk of infection**

- ▶ Subject the device to a hygiene treatment at regular intervals as described in this section.

### **NOTICE!**

#### **Damage to device**

Unsuitable cleaning agents and disinfectants may damage the sensitive surfaces of the device and impair its operability.

- ▶ Do not use aggressive or abrasive cleaning agents.
- ▶ Do not use organic solvents (e.g. white spirit or petroleum spirit).

### 7.1 Cleaning

- ▶ If required, moisten a soft cloth with a mild soap solution and wipe the device over with it.

### 7.2 Disinfecting

1. Disinfect the device at regular intervals using a disinfectant suitable for sensitive surfaces and acrylic glass. seca recommends a 70 % ethanol solution as a surface disinfectant.
2. Follow the instructions for use of the disinfectant.
3. Disinfect the device:
  - ▶ Moisten a soft cloth with disinfectant and wipe down the device with it.
  - ▶ Comply with the intervals, see table:

| Interval                                        | Component                                                                                                    |
|-------------------------------------------------|--------------------------------------------------------------------------------------------------------------|
| <b>Before</b> and <b>after</b> each measurement | Tray                                                                                                         |
| As required                                     | <ul style="list-style-type: none"> <li>• Display</li> <li>• Controls</li> <li>• Weighing platform</li> </ul> |

### 7.3 Sterilizing

The device must not be sterilized.

## 8 FUNCTION CHECK

### CAUTION! Personal injury

If you find faults or deviations during the function check which you are unable to remedy with the aid of the “Troubleshooting” section in this document, you must not use the device.

- Have the device repaired by seca Service or an authorized service partner.
- Follow the section entitled “Servicing” in this document.

- Perform a function check before each use.

A complete function check includes:

- A visual inspection for mechanical damage
- A test of the alignment of the device
- A visual and function check of the display elements
- A function check of all the controls shown in the section entitled “Overview”
- A function check of the optional accessories

If you find faults or deviations during the function check, initially try to remedy the fault with the aid of the “Troubleshooting” section in this document.

## 9 TROUBLESHOOTING

### 9.1 General faults

| Fault                                                                                            | Cause                                                      | Remedy                                                                                                                                                                                                                                                                                                        |
|--------------------------------------------------------------------------------------------------|------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| No weight is displayed even though there is a load on the scale                                  | Device has no power supply                                 | <ul style="list-style-type: none"><li>• Switching on the device</li><li>• In battery mode: Inserting batteries</li><li>• In power supply operation: Connecting the power supply unit</li></ul>                                                                                                                |
| <b>0.000</b> is not displayed before weighing takes place                                        | A load was already on the device before it was switched on | <ul style="list-style-type: none"><li>• Remove the load from the device</li><li>• Switch off the device, then switch it back on again</li></ul>                                                                                                                                                               |
| Weight display flashes                                                                           | Consistent weight value has not yet been detected          | <ul style="list-style-type: none"><li>• One of the following functions is active: <b>Hold</b>, <b>Tare</b>. Wait until the device has detected a consistent weight value. The display stops flashing.</li><li>• Take the load off the device, wait until <b>0.000</b> is displayed and weigh again.</li></ul> |
| One display segment is displayed continuously or not at all                                      | Display is defective                                       | <ul style="list-style-type: none"><li>• Inform seca Service</li></ul>                                                                                                                                                                                                                                         |
|  is displayed | Battery voltage is dropping                                | <ul style="list-style-type: none"><li>• Replace the batteries</li></ul>                                                                                                                                                                                                                                       |
| <b>bAtt</b> is displayed                                                                         | Batteries have discharged                                  | <ul style="list-style-type: none"><li>• Replace the batteries</li></ul>                                                                                                                                                                                                                                       |
| <b>StOP</b> is displayed                                                                         | Maximum capacity has been exceeded                         | <ul style="list-style-type: none"><li>• Remove the load from the device</li></ul>                                                                                                                                                                                                                             |
| - - - - is displayed before weighing takes place                                                 | A load was already on the device before it was switched on | <ul style="list-style-type: none"><li>• Remove the load from the device</li><li>• Switch off the device, then switch it back on again</li></ul>                                                                                                                                                               |

17-10-07-688-100a\_2026-09S

| Fault                           | Cause                                                                      | Remedy                                                                                                                                                                                                                                    |
|---------------------------------|----------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Er:[No.]:11</b> is displayed | The device has too high a load or too high a load in one corner            | <ul style="list-style-type: none"> <li>Remove the load from the device or distribute the weight more evenly</li> <li>Switch off the device, then switch it back on again</li> <li>If the problem persists, inform seca Service</li> </ul> |
| <b>Er:[No.]:12</b> is displayed | The device has been switched on with too high a load                       | <ul style="list-style-type: none"> <li>Remove the load from the device</li> <li>Switch off the device, then switch it back on again</li> </ul>                                                                                            |
| <b>Er:[No.]:16</b> is displayed | The device was caused to oscillate, the zero point could not be determined | <ul style="list-style-type: none"> <li>Switch off the device, then switch it back on again</li> </ul>                                                                                                                                     |

## 10 SERVICING

### CAUTION!

**Faulty measurements as a result of verification being performed poorly or not at all**

- Have verification performed only by authorized persons.
- Always have verification performed if one or more quality seals are damaged.

seca recommends having your device serviced prior to verification.

### CAUTION!

**Faulty measurements as a result of poor servicing**

- Have servicing and repairs carried out exclusively by seca Service or an authorized service partner.
- You can find your local service partner at [www.seca.com](http://www.seca.com).

## 11 TECHNICAL DATA

### 11.1 General technical data

| General technical data        |                                     |
|-------------------------------|-------------------------------------|
| Dimensions                    |                                     |
| • Depth                       | 362 mm                              |
| • Width                       | 550 mm                              |
| • Height                      | 165 mm                              |
| Net weight                    | approx. 2.5 kg                      |
| Ambient conditions, operation |                                     |
| • Temperature                 | +10 °C to +40 °C / +50 °F to 104 °F |
| • Air pressure                | 700 – 1060 hPa                      |
| • Humidity                    | 30 % – 80 % no condensation         |
| Ambient conditions, storage   |                                     |
| • Temperature                 | -10 °C to +65 °C / +14 °F to 149 °F |
| • Air pressure                | 700 – 1060 hPa                      |
| • Humidity                    | 0 % – 95 % no condensation          |

| General technical data                                                                                                                                                 |                                                                                     |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|
| Ambient conditions, transport <ul style="list-style-type: none"> <li>• Temperature</li> <li>• Air pressure</li> <li>• Humidity</li> </ul>                              | -10 °C to +65 °C / +14 °F to 149 °F<br>700 – 1060 hPa<br>0 % – 95 % no condensation |
| Height of digits                                                                                                                                                       | 22 mm                                                                               |
| Power supply <ul style="list-style-type: none"> <li>• Batteries <ul style="list-style-type: none"> <li>– Supply voltage</li> <li>– Battery type</li> </ul> </li> </ul> | 9 V<br>6 x type AA, 1.5 V                                                           |
| Current consumption                                                                                                                                                    | Approx. 30 mA                                                                       |
| Maximum runtime in battery mode                                                                                                                                        | Approx. 60 hours                                                                    |
| Medical device in accordance with regulation (EU) 2017/745                                                                                                             | Class I with measuring function                                                     |
| EN 60601-1: Medical electrical equipment, type B                                                                                                                       |  |
| Type of protection                                                                                                                                                     | IP20                                                                                |
| Duty cycle                                                                                                                                                             | Continuous duty                                                                     |
| Application parts in accordance with IEC 60601-1                                                                                                                       | Tray, weighing platform                                                             |

## 11.2 Weight measurement

| seca 384                                                                                                                                                                                                                                              |                                     |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------|
| Verification in compliance with directive 2014/31/EU                                                                                                                                                                                                  | Class III                           |
| Maximum capacity <ul style="list-style-type: none"> <li>• Weighing range 1</li> <li>• Weighing range 2</li> </ul>                                                                                                                                     | 10 kg<br>20 kg                      |
| Minimum capacity <ul style="list-style-type: none"> <li>• Weighing range 1</li> <li>• Weighing range 2</li> </ul>                                                                                                                                     | 0.2 kg<br>0.4 kg                    |
| Graduation <ul style="list-style-type: none"> <li>• Weighing range 1</li> <li>• Weighing range 2</li> </ul>                                                                                                                                           | 10 g<br>20 g                        |
| Tare range <ul style="list-style-type: none"> <li>• Weighing range 1</li> <li>• Weighing range 2</li> </ul>                                                                                                                                           | 10 kg<br>20 kg                      |
| Accuracy on initial verification <ul style="list-style-type: none"> <li>• Weighing range 1: 0.2 kg – 5 kg</li> <li>• Weighing range 1: 5 kg – 10 kg</li> <li>• Weighing range 2: 0.4 kg – 10 kg</li> <li>• Weighing range 2: 10 kg – 20 kg</li> </ul> | ± 5 g<br>± 10 g<br>± 10 g<br>± 20 g |

| seca 385                                                                                                          |                |
|-------------------------------------------------------------------------------------------------------------------|----------------|
| Verification in compliance with directive 2014/31/EU                                                              | Class III      |
| Maximum capacity <ul style="list-style-type: none"> <li>• Weighing range 1</li> <li>• Weighing range 2</li> </ul> | 20 kg<br>50 kg |
| Minimum capacity <ul style="list-style-type: none"> <li>• Weighing range 1</li> <li>• Weighing range 2</li> </ul> | 0.4 kg<br>1 kg |

| seca 385                           |        |
|------------------------------------|--------|
| Graduation                         |        |
| • Weighing range 1                 | 20 g   |
| • Weighing range 2                 | 50 g   |
| Tare range                         |        |
| • Weighing range 1                 | 20 kg  |
| • Weighing range 2                 | 50 kg  |
| Accuracy on initial verification   |        |
| • Weighing range 1: 0.4 kg – 10 kg | ± 10 g |
| • Weighing range 1: 10 kg – 20 kg  | ± 20 g |
| • Weighing range 2: 1 kg – 25 kg   | ± 25 g |
| • Weighing range 2: 25 kg – 50 kg  | ± 50 g |

## 12 OPTIONAL ACCESSORIES AND SPARE PARTS

| Accessory/spare part                        | Article number |
|---------------------------------------------|----------------|
| Carry case <b>seca 413</b>                  | 413 0000 009   |
| Carry case <b>seca 414</b>                  | 414 0000 009   |
| Measuring board <b>seca 417</b>             | 417 1721 009   |
| Carry case <b>seca 431</b>                  | 413 0000 009   |
| Baby scale cart <b>seca 402</b>             | 402 0000 009   |
| Baby scale cart with drawer <b>seca 403</b> | 403 0000 009   |

## 13 DISPOSAL

### 13.1 Disposing of the device

Do not dispose of the device in your household waste. The device must be properly disposed of as electronic scrap. Follow your respective national regulations. For more information, please contact seca Service or an authorized service partner.

### 13.2 Disposing of batteries and rechargeable batteries

Do not dispose of used batteries and rechargeable batteries in household waste, regardless of whether they contain harmful substances or not. As a consumer, you are legally obliged to dispose of batteries and rechargeable batteries via local authority collection points or trade collection points. Only dispose of batteries and rechargeable batteries once they are completely discharged.

## 14 WARRANTY

There is a two-year warranty period from delivery for defects attributable to poor materials or workmanship. All movable parts, e.g. batteries, cables, power supply units, rechargeable batteries etc. are exempt. Defects which come under the warranty will be repaired for the customer free of charge against proof of purchase. Additional claims cannot be considered. Costs of transport to and from seca are the responsibility of the customer if the device is located somewhere other than the customer's headquarters. In the event of transport damage, claims under warranty can only be made if the complete original packaging was used for transport and the device was secured and fastened in it according to its originally packaged condition. You should therefore keep all packaging parts.

The warranty will be voided if the device is opened by persons not expressly authorized by seca to do so.

In the event of a warranty issue, please contact your local seca office or the dealer from whom you ordered the product.

## 15 DECLARATION OF CONFORMITY

seca gmbh & co. kg hereby declares that the product complies with the terms of the applicable European directives and regulations. The unabridged declaration of conformity can be found at [www.seca.com](http://www.seca.com).

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